

ORCHESTRA Cool Cryolipolysis Syetem

Fat treatment & Cryotherapy

Abdomen

Thight

Chest

Love Handle

Buttock

Upper Back

Arm


Leg



Working priciple


As triglyceride in fats will be converted into solid in particular low temperature, this machine uses precisely controlled cooling energy to the target areas for treatment. It selectively targets at certain fat bulges and eliminate fat cells through a gradual process that does not harm the surrounding tissues and reduce unwanted fat, When fat cells are exposed to precisely control cooling, the triglyceride in fats will be converted into solid, and they will be Crystallization, aging, and then die. The dead fat cells will be out of body by natural metabolism.

Treatment Process

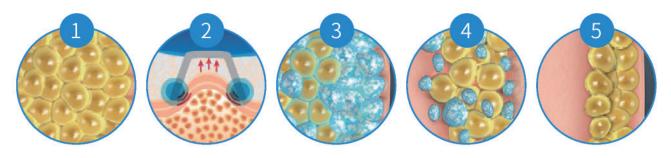

- 01 | Treatment area prior to applying cooling
- 02 | Precisely controlled cooling is applied to the treatment area
- 03 Underlying fat cells are targeted
- 04 | Fat cells begin a natural removal process that continues for several months
- 05 l The natural removal of fat cells over time results in gradual fat reduction

Dynamic handpiece

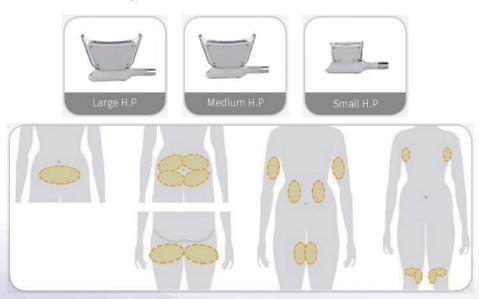

- Four different sizes handpieces optional, easy treatment for whole body.
- Four handpieces working simultaneously, save staff cost, save treatment time.
- → Unique double chin treatment handpiece, lead the latest technology in the world.
- ⇒ Super strong cooling system ensure the minimun treatment temperature reduce to -15 → 5°C
- Non-surgery, no risk, easily intelligentize and humanization operation.
- No need recovery, no downtime, come back to work immediately.
- > Variety accessories for multi-functional treatment.
- Luxury design.
- → Easy operation for all program.



Specifications

Handles number	4 handles (1 Large, 2 Medium, 1 Small)
Handle screen	5 inch touch screen
Vacuum intensity	0-100kPa adjustable
Handle sizes	Large : 220*76*125mm; Medium :160*56*65mm; Small : 125*45*70mm
Machine screen	10.4 inch touch screen
Cooling temperature	-15℃5℃
Cooling system	Water, Air, Semi-conductor triple cooling system
Input power	2500W
Output power	1600W
Machine Dimensions	1230*450*420mm
Fuse specifications	Ø5×25 15A
Machine Dimensions	AC220V±10%,10A;50HZ/ AC110V±10%10A,60HZ





world gate medical

Erbil:

0751 235 4023

Dohuk:

0750 303 2672

Kirkuk:

0750 204 7132

Baghdad:

0750 835 0758

Diyala:

0772 511 1466

Sulaimanya:

0750 852 8458

Manager:

0751 441 1855



Service: 0751 434 3415 & 0782 838 6448

www.worldgate-medical.com sales@worldgate-medical.com
Main Office: Iraq.Erbil.Mahmood Alhafeed Square.Qessab building